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TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY. —

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eta 0706

INSTRUCTIONS FOR USE

Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the warranty, the receipt and, if possible, the box with the internal packing.

I. GENERAL SAFETY INSTRUCTIONS

- Consider the instructions for use as a part of the appliance and pass them on to any other user of the appliance.
- Please check if the voltage specified on the rating plate complies with the voltage in your electrical socket (this product is designed to operate within a range of 100-240 volts). The charging base must be connected to electric socket in the electrical mains complying with the applicable standards.
- Children and persons with limited physical or mental capabilities or lack of experience may only use this appliance if they are supervised or have been instructed of the use of the appliance in a safe manner and understand possible risks. Unattended children must not perform cleaning and maintenance intended to be performed by users. Children must not play with the appliance.
- Before replacing accessories or accessible parts, which move during operation, before assembly and disassembly, before cleaning or maintenance, turn off the appliance and disconnect it from the mains by pulling the charger from the power socket!
- Never use the appliance if the USB power cable or adapter is damaged, if it does not work properly, if it fell to the ground and got damaged. In this case take the appliance to a special service to check its safety and proper function.
- The appliance has no charging source of its own, it is only designed for charging via USB (an adapter with USB output, with output voltage of 5.0 V, and output current of at least 700 mA). Before making any connection, please pay attention to instructions of the manufacturer of the respective power source.
- The trimmer with the plugged adapter must not be used at places where it could fall into a bath tub, a wash basin or a swimming pool. If the trimmer fell to water anyway, do not take it out! First of all, unplug the adapter from the socket and only then remove the trimmer. In this case take the appliance to a special service to check its safety and proper function.
- Do not leave the iron plugged in power supply without supervision!
- Always unplug the appliance from charger if you leave it unattended or before cleaning.

- The appliance is only intended for individual use at home for personal oral cavity care! It is not intended for commercial use, for more than one patient in dentist's consulting rooms or in any other healthcare facilities!
- **Do not put or take out the plug to or from a socket with wet hands!**
- Use only the accessories designed for this type. Do not use damaged adapters.
- Do not use the application in exterior or where aerosol products are used (**such as sprays**), or where oxygen is handled.
- Guarantee of 24 months does not apply to decrease in the accumulator's effective capacity caused by its use or age. Effective capacity gradually decreases depending on the method of use of the accumulator.
- The appliance is not designed to care for animal teeth.
- Only use the charging base at horizontal position at places where overturning is not imminent, and in a sufficient distance from heat sources (such as **stove, fireplace, or heat radiators**).
- Clean the accessories after each use. Failure to adhere to this instruction might cause worsening your state of health, because without disinfection, spread of mildew and bacteria is supported.
- Do not insert the appliance or its accessories in any body cavities, except for oral cavity.
- Charge the accumulator at the normal room temperature. Do not expose the accumulator to temperatures higher than 50 °C to prevent causing damage to the accumulator.
- The accumulator must not be removed when the appliance is normally used. Only remove the accumulator if the appliance is to be discarded (see article **V. ECOLOGY**).
- **Do not leave the appliance unattended!**
- Protect the appliance from dust, chemicals, toiletries, liquid cosmetics, large fluctuations in temperature, and direct sunshine.y.
- Do not attempt to repair the appliance by yourselves.

SPECIAL SAFETY INSTRUCTIONS

- If you have unpleasant feeling in your mouth or if you feel pain during the cleaning, stop using the appliance and contact your dentist. Contact the dentist also in the case that you have any doubt about use of the appliance in your individual case.
- If your toothpaste contains peroxide, baking soda or other bicarbonate (frequently used in whitening toothpastes), thoroughly clean the brush and its handle after use to prevent damaging the appliance.
- Do not use the appliance in care for children, incapacitated persons or persons who are unable to express their feelings.
- **If you have undergone tooth or gum operation recently, suffered long-term problems with your teeth or were subject to extensive medical interventions in the last 2 months, or if you suspect problems of dental character (with fillings, crowns, etc.), consult your dentist before using the appliance.**
- **If you experience any health troubles (for example considerable bleeding, or if any bleeding persists even after 2 weeks of use of the brush), stop using the appliance and contact your physician for consultation.**
- For hygienic reasons do not share the toothbrush with another person or more persons. Infections could spread or oral cavity inflammation could be initiated.
- The appliance complies with the safety standards for electromagnetic radiation. However if you use any of the following devices: cardiostimulator, electrocardiograph or other medical device implanted in your body or used as a support, consult your physician or the implant manufacturer before use.
- The power cord of the adapter must not be damaged with sharp or hot objects, open fire and it must not be sunk into water.

- The cable supplied with this appliance is intended for this appliance only, do not use it for other purposes. Also, use only the supplied cable with this appliance, do not use any other cable for charging (e.g. from another appliance).
- Only use the appliance with the original accessories designated for the respective type. Do not use damaged/deformed adapters.
- Only use this appliance including its accessories for the designated purpose as described in this instruction manual. The appliance is only designed for cleaning teeth, gums and tongue. Never use the appliance for any other purpose.
- **WARNING:** There is a risk of injury in case of incorrect use (not in accordance with the manual) of the appliance.
- All text in other languages and pictures shown on the packaging, or the product itself, are translated and explained at the end of this language mutation.
- The producer does not undertake any responsibility for damage caused by improper use of the appliance (e.g. **injury etc.**) and its guarantee does not cover the appliance in the case of failure to comply with the safety warning above.

II. DESCRIPTION OF FEATURES (FIG. 1)

A – Brush body

- A1 – ON/OFF switch  + mode switch
- A2 – battery/active mode indicator light
- A3 – protective cover
- A4 – Power USB slot

B – Brush (2 pcs included)

- B1 – brush head
- B2 – coloured ring to distinguish users

C – USB cable * (Use AC power adapter with USB or mobile power)

Note: Replacement brushes **B** can be purchased under order number ETA 070690100.

III. INSTRUCTIONS FOR USE

Remove all packaging materials and take out the appliance and any accessories. Parts which come in contact with the oral cavity must be rinsed in warm drinking water before first use.

CHARGING THE BRUSH (Fig. 2).

- 1) Make sure the appliance is off. While charging, the appliance cannot be in use. If the appliance is turned on and you start charging, the appliance shall switch off.
- 2) Put the provided charge USB cable **C** into the charge slot **A4** of brush body, connect the other charger compatible AC wall adapter or mobile power.
- 3) The charging procedure is signalled by a flashing **A2** red LED indicator.
- 4) When the battery is fully charged, the **A2** LED indicator it will go off. Unplug the adaptor from the wall socket first and only then unplug the connector from the case.
- 5) For a full charge, the battery needs to be charged for approximately 4 hours. Depending on the cleaning mode used, a fully charged battery will provide approximately 80 minutes of operation.

Note: Although the brush should not be charged before it is fully discharged, we recommend charging it at regular intervals (e.g. after 2-3 weeks of use) in order to ensure optimum battery life.

Caution

- During operation, the **A2** LED indicator is lit, when the battery of the accumulator is low and it needs to be recharged, the LED indicator starts flashing (10x).
- **Do not exceed the recommended time of charging during regular charging!**
- The appliance warms up during charging, but this is normal and not a defect.
- The operating time is only approximate. The time depends on the condition of the battery, ambient temperature, age, method of charging, etc.

USE

Before use, slide the brush **B** onto the brush body **A** to connect them (See Fig. 1). A small gap will remain between the brush body **A** and the coloured ring **B2** which is not a defect. To remove the brush, proceed the other way round. Do not switch on the appliance without the brush attached!

Moisten the brush and apply a small amount of toothpaste on it. Put the brush to the teeth at a slight angle towards the gums (see Fig. 3) and switch on any mode firstly by pressing the on/switch button **A1**. Move the brush slowly along the teeth back and forth, applying a slight pressure, so as not to miss any tooth or area in the mouth. Clean each tooth for about 2-3 seconds.

The brush operating time is limited to 2 minutes which corresponds to the time for cleaning recommended by experts. If you wish to continue cleaning, simply switch on the brush again. Each part of the mouth (quadrant) should be cleaned for about 30 seconds (see Fig. 4). A short interruption of the brush operation will prompt you to move to another quadrant after this time elapses.

Notes

- If this is your first sonic toothbrush and you have not had any experience with this cleaning technology, it might take you a while to get used to the specific vibration (tingling). However, most users say that they got used to it quickly after a few uses and cleaning became perfectly natural for them.
- Do not apply excessive pressure on your teeth.
- You can also use this brush when you have fixed braces, facets, crowns, etc. However, please bear in mind that if the brush is used with braces, the brush bristles are subject to more wear.
- We recommend that you replace the brush **B** with a new one after 3 months of use. However, if the brush is showing signs of significant wear or damage.

Mode on/switch button control **A1**

When you press the button once, the brush switches on in mode **C**. If you press the button again within 5 seconds, the brush will switch to mode **S**. If you press the button repeatedly within 5 seconds, the brush will switch to mode **W**. If you press the button repeatedly within 5 seconds, the brush will switch to mode **M**. When you press the button again, the brush will switch off. The brush will also switch off if you press the button **A1** in mode **C/S/W** after 5 seconds after its activation.

Cleaning modes:

- 1) "**C**" – **CLEAN** mode (blue light **A2** is on)
 - standard intensive mode for regular cleaning of healthy teeth.

- 2) "S" – **SOFT** mode (green light **A2** is on)
 – gentle mode for sensitive teeth or mouth parts; recommended in case of increased sensitivity of teeth.
- 3) "W" – **WHITEN** mode (red light **A2** is on)
 – cleaning mode adapted for better teeth whitening,
- 4) "M" – **MASSAGE** mode (blue light **A2** is flashing)
 – pulse mode for gentle gum massaging.

Activation of the mode selected is indicated by the indicator **A2** coming on.

IV. CLEANING AND MAINTENANCE

For maximum hygiene, clean the brush B after each use! Rinse the brush head **B1** and the brush body **A** under running water and wipe the brush body with a damp cloth. Use the damp cloth to wipe the charging base and its cable as well. Never immerse the brush body **A** in water! Do not use a dishwasher for cleaning! Allow the brush body **A** to dry properly before re-use. Store the appliance and all the accessories in a dry, dust-free and safe place, out of reach of children and incapacitated persons.

V. ENVIRONMENTAL PROTECTION

If the dimensions allow, there are labels of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal of the product, hand it over at special collection places where it will be taken over free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal. Fines can be imposed for improper disposal of this type of waste in agreement with the national regulations. This appliance is equipped with Li-ion batteries with a long service life. For environmental protection it is necessary to remove the battery from the appliance after its service life and the battery and the appliance must be disposed of safely using special collection networks. Ask for more details at the local authority or at the local waste collection site.

Removing the battery

Remove the battery from the appliance only when completely discharged (i.e. leave the appliance running until it stops).

The procedure for removing the battery is shown in Fig. 5:

1. Using a screw-driver, remove the screw cover and loosen the screw.
2. Hold the brush body **A** upside down, press the charging part (see arrows in the figure) and remove the 2 hooks which hold the inner components of the brush body (**A**).
3. Remove the battery from the compartment using a screw-driver.

VI. TECHNICAL DATA

Voltage (V)	shown on the type label of the appliance
Weight (kg) approximately	0.11
Size, (mm):	180 x 27 x 27
Noise level: Acoustic noise level of ≤ 65 dB(A) re 1pW.	

The manufacturer reserves the right to change the technical specifications and accessories for the respective models. Expert service must do the large maintenance or to intervene into inside parts of the appliance! Non-observance of the instructions of the manufacturer cancels the right for guarantee repair!


*HOUSEHOLD USE ONLY; DO NOT IMMERSE IN WATER OR OTHER LIQUIDS;
TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY
FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS,
BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.*

NOTICE

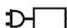- Please read the instructions manual

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© DATE 06/05/2021