<i>IPL epilátor</i> • NÁVOD K OBSLUZE	CZ	3-16
<i>IPL epilátor</i> • NÁVOD NA OBSLUHU	SK	17-30
<i>IPL epilator</i> • INSTRUCTIONS FOR USE	GB	31-43
<i>IPL epilátor</i> • KEZELÉSI ÚTMUTATÓ	H	44-57
<i>IPL epilator</i> • INSTRUKCJA OBSŁUGI	PL	58-71
<i>IPL Epilierer</i> • BEDIENUNGSANLEITUNG	D	72-86




Fenité

18/2/2020

eta

1

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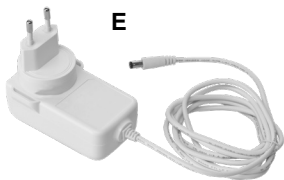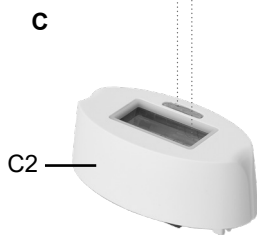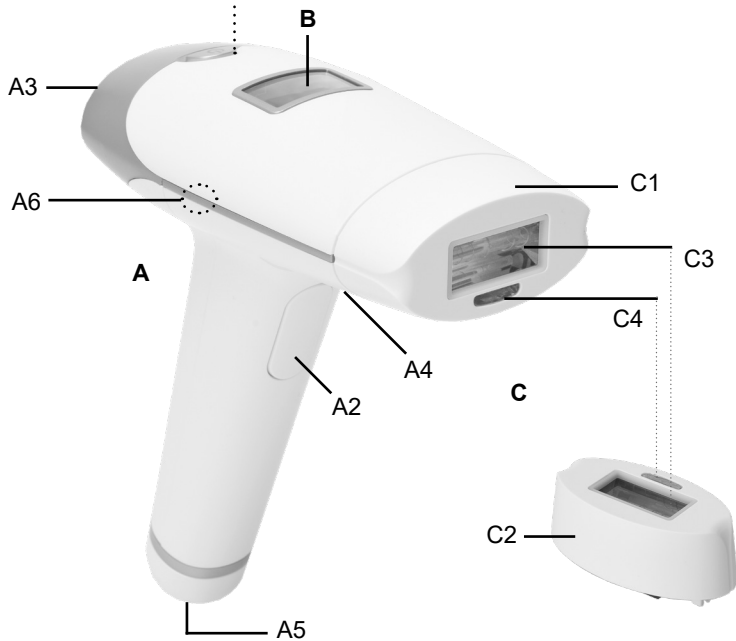TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY. *Nebezpečenstvo udusenía. Nepoužívajte toto vrecúško*

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UPOZORNENIE

INSTRUCTIONS FOR USE

Before using this device, please read carefully following instructions, including all contraindications, warnings and safety information and keep them for later use.

Due to its specific technology, the device is not suitable for everyone. The most effectively it works on light to medium skin tones and dark hair. Therefore, please get familiarized with following information, including the skin tone and hair colour tables, to make sure that its use is suitable and safe for your skin.

I. SAFETY WARNING

- Before the first use, read the instructions for use carefully, look at the picture and keep the instructions for use.
- Check whether the data on the type label correspond with the voltage in your socket.
- **This appliance must not be used by children. Keep the device and its supply cord out of reach of children. The device may be used by persons with reduced physical or mental abilities or lack of experience and knowledge only if they are under supervision or they have been instructed about using the appliance safely and understand the potential dangers. Children must not play with the appliance.**
- Before replacing accessories or accessible parts, before assembly and disassembly, before cleaning or maintenance, turn off the appliance and disconnect it from the mains by pulling the charger from the power socket!
- **Do not leave the iron plugged in power supply without supervision!**
- If the power cord of the charger is corrupted, it has to be replaced by the manufacturer, its service technician or a similarly qualified person so as to prevent dangerous situations.
- Never operate this device if its adapter, cord, plug, attachments or their glass filters are damaged, if it does not work properly, if it has been dropped and damaged, or dropped into water. In such cases bring the device to an authorized service centre to have its safety and proper operation examined.
- The appliance with the plugged adapter must not be used at places where it could fall into a bath tub, a wash basin or a swimming pool. If the appliance fell to water anyway, do not take it out! First of all, unplug the adapter from the socket and only then remove the appliance. In this case take the appliance to a special service to check its safety and proper function.
- Always unplug the appliance from charger if you leave it unattended.

- For additional protection we recommend installation of current protection (RCD) with the nominal release current not exceeding **30 mA**. Ask a control technician or an electrician for advice.
- **WARNING** – Certain parts of the device (glass filters of the attachments) may become very hot and cause burns. Do not touch them during operation. Special attention should be paid to the presence of children or persons with reduced physical, sensory or mental capabilities.
- **Considering hormonal and other changes in puberty, the device is not intended for use by persons under 18 years of age.**
- **Facial attachment C2 is suitable for persons above 25 years of age. However, younger persons, due to growth dynamics and hormonal changes, are not recommended to use this attachment!**
- Never switch on the device without the attachment placed. Do not use the attachment if the glass filter, which forms its integral part, is missing. At the same time, check regularly, if the glass filter is not damaged.
- Do not use the device for any other surface than human skin (excluding forbidden areas). Also do not use it, if the application window is not in contact with your skin.
- Leave the ventilation slots free and in no way block or cover them during the operation. Ensure good circulation of air in the room.
- Use the device only with the intensity setting, which is suitable for your skin tone/type (see below). If a higher intensity than that allowed is used, the risk of side effects and undesired reactions of the skin may increase.
- To prevent the risk of skin injury, do not treat the same area with flashes (nor if changing the attachments!) more than once a week. Do not treat the same area with more flashes during one session. Repeated treatment due to different hair growth stages does not improve the effectiveness of epilation (as well as the treatment with facial attachment) but occurs the risk of skin damage.
- Never look directly at the flash produced by the device (at the glass filter); nor try to switch on the device when pointing to eyes, even with safety eyewear on.
- Never use the device, if small children or toddlers are present in the room! The light produced may be unpleasant to them.
- Before epilation, it is necessary to remove all metal objects, such as rings, bracelets, necklets, piercings etc. You must not use any cosmetic products or, for example, make-up pencils containing metal to mark the places to be treated! Do not use it, for example, for lips contoured with a permanent or semi-permanent contour pencil.
- Before use, never clean the areas to be treated with agents containing alcohol (alcohol, acetone etc.).
- In any case of activation (flashing) of the device, even if the touch sensor with application window is not in contact with the skin, do not use the device and have it examined by an authorized service centre.
- The device is intended for use at home and for similar purposes! It is not designed for use in cosmetic salons or for any other commercial purposes!
- **Do not put or take out the plug to or from a socket with wet hands!**
- **Never immerse the device, its attachments or adapter in water or other liquids (nor partially), do not expose them to moist and do not wash them under the tap!**
- Use only the accessories designed for this type. Do not use damaged adapters.

- The appliance must not be used in damp or wet environment and in any environment with the danger of fire or explosion (spaces where chemicals, fuels, oils, gases, paints and other flammable or volatile materials are stored).
- The power cord of the adapter must not be damaged with sharp or hot objects, open fire and it must not be sunk into water.
- The adaptor supplied with this appliance is intended for this appliance only, do not use it for other purposes. Also, use only the supplied adaptor with this appliance, do not use any other adaptor for using (e.g. from another appliance).
- Never use the appliance for any other purpose than for the intended purpose described in these instructions for use!
- The producer is not responsible for any loss caused by improper handling of the appliance and its accessories (e. g. **injury, appliance damage, fire**, etc.) and its guarantee for the appliance does not apply in the case of non-compliance with the safety instructions above.

II. CONTRAINDICATIONS AND SPECIFICS OF USE

Contraindications

This type of device is unique with the IPL technology used. It allows to bring the high professional effectiveness to the comfort of your home. However, this technology has some limitations and specifics of use – it is not intended for all groups of persons and it cannot be used under certain circumstances or in case of certain diseases.

VERY IMPORTANT

Any failure to follow the instructions and recommendations given here may result in skin damage, aggravation of possible current problems, or occurring of adverse side effects, such as burns, blisters, discolorations or scarification of the skin.

If you have any doubts concerning the use of the device in your individual case consult your dermatologist or general practitioner first! Further, in case of your long-term (in certain cases also short-term) use of certain medication (or hormonal contraception), ask if they may not influence the sensitivity of your skin to light and heat and thus negatively influence the skin reaction to IPL treatment!

The device **must not** be used by:

- Children and elder persons;
- Persons with local or systemic inflammatory disease;
- Persons with implanted medical devices and systems (such as pacemaker, skin ports etc.).
- Persons with precancerous lesions or several atypical moles on the area to be treated with the device.
- Persons with skin diseases (e.g. skin cancer, psoriasis, vitiligo), with medical history of skin cancer or any other focal lesions of cancer localised on the areas of intended use of the device or the persons currently being treated with chemotherapy or radiotherapy.
- Persons with photosensitive skin or taking photosensitizing medication;
- Persons with any chronic disease, heart disease, thyroid gland disease, further the persons with the diseases such as lupus erythematosus, diabetes, epilepsy or retina disease or receiving an immunosuppressive therapy.
- Pregnant or breastfeeding women or women attempting to conceive.
The device has not been tested for this group of persons.

Further, do not use the device in following cases:

- If the skin on the area to be treated has been treated with peeling, sunburned, cut, has open wounds, eczemas, lesions, haematomas, burns, inflammation of hair follicles, abrasions or is currently being treated with injections or if you have had a surgery or facial rejuvenation surgery in given area (such as Botox, facial filler, thread lifts etc.).
- If you have varicose veins, dilated veins or a history of any other vascular disorder on the area to be treated.
- If you have a history of collagen disorder, including a history of keloid scar formation or a history of poor wound healing, or you have skin implants.
- If you have any bleeding disorder; also, if you use anticoagulation medications.
- Over or near silicone implants or piercings.
- Immediately after cosmetic treatment (e.g. after AHA, botulotoxin treatment etc.).
- If you visit tanning centre.
- In moles (neither near them nor over them!), freckles, darker pigmented areas, scars or other skin irregularities.
- If you have warts, freckles, tattoos or permanent make-up on the area to be treated
- If you have used long-lasting deodorant on the area to be epilated.
- If your skin tone is very dark (see part “Skin type „below). The dark skin absorbs more light energy, which may result in unpleasant or painful sensation and undesired effects (e.g. burning, blisters formation, discoloration or scarification) and damage to skin.
- If your skin is very sensitive.
- If you suffer from solar urticaria (allergy to sun).
- Before or shortly after IPL application, do not use any cosmetic products (e.g. products with pigments, perfumes, alcohol etc.).
- In any other case which would be considered by your doctor as risky for the use of this device.

The above-mentioned lists are not complete and do not include all possible cases and complications, that may exclude the use of IPL technology. In case of any doubts you should contact your dermatologist or general practitioner.

Further, special attention should be paid to the issue of tanning. Do not use the device on recently tanned or artificially tanned skin! All types of tanning (natural, as well as artificial one, including the use of self-tanners) may reduce the safety of IPL treatment, as the skin becomes very sensitive. Do not use also if you are planning to travel soon e.g. to the sea and to get exposed to the sunlight. It applies also vice versa - after IPL treatment, do not expose your unprotected skin to any tanning or sunlight for at least 2-3 weeks. Regarding the above-mentioned, it is not recommended to use the device in summer months (or in days with increased sunlight, regardless the season of the year), when people are generally very often exposed to the sunlight.

Skin type

IPL technology has its advantages and limitations. One of the main factors, determining the resulting effect and the safety of use of the IPL device, is the skin tone and the colour of hair. To put it simply, the best results are achieved by the persons with light skin and dark hair. The darker the hair, the more energy from the flashes is absorbed, making the epilation more effective. However, the darker skin tone is less desirable from the point of view of the safety in case of skin discoloration.

Therefore, the device is not intended for use by the persons with very dark skin tone, as there is a risk of skin injury (e.g. burning). If your skin tone is very dark, do not use the device.

If you have a darker skin tone (and before using the device for the first time) always start the epilation at the lowest intensity. Only after you make sure that given setting does not cause any complications, you can increase the light intensity setting, if needed.

Please also note that IPL device is not effective in persons with light blond, white, red or grey hair, since due to the low level of dark pigment in this case, the effectiveness of the epilation may be very low.

Hair colour	Skin tone type					
	1	2	3	4	5	6
Black	✓ ■ ■ ■ ■ □	✓ ■ ■ ■ □ ⊗	✓ ■ ■ □ □ ⊗	✓ ■ □ □ ⊗ ⊗	✓ ■ □ ⊗ ⊗ ⊗	✗
Dark brown	✓ ■ ■ ■ □ □	✓ ■ ■ ■ □ ⊗	✓ ■ ■ □ □ ⊗	✓ ■ □ □ ⊗ ⊗	✓ ■ □ ⊗ ⊗ ⊗	✗
Brown	✓ ■ ■ ■ □ □	✓ ■ ■ ■ □ ⊗	✓ ■ ■ □ □ ⊗	✓ ■ □ □ ⊗ ⊗	✓ ■ □ ⊗ ⊗ ⊗	✗
Dark blond /light brown	✓ ■ ■ ■ □ □	✓ ■ ■ ■ □ ⊗	✓ ■ ■ □ □ ⊗	✓ ■ □ □ ⊗ ⊗	✓ ■ □ ⊗ ⊗ ⊗	✗
Light blond	✗	✗	✗	✗	✗	✗
Red	✗	✗	✗	✗	✗	✗
White / grey	✗	✗	✗	✗	✗	✗

- ✓ - can be used; ✗ - do not use
- - recommended intensity
- - the intensity can be used
- - the intensity can be used with high precaution
- ⊗ - the intensity cannot be used

The above displayed table shows the manner of use of the device and recommended setting for given skin/hair type. To put it simply, the darker skin you have, the lower light intensity should be set. If your skin type is approximately 6 or darker, do not use the device. Also do not use the device in case you do not have standard skin pigmentation (e.g. vitiligo or albinism). Determine your skin tone according to the reference zone (skin area), which has not been exposed to the sunlight – ideally the area hidden under bikini or underarms. The colour tones including the table are illustrative only and serve only for better understanding of the device operation in relation to various skin and hair colour tones. Always respect the warnings in the instructions, including the skin sensitivity test.

Never use the device on following body or facial areas:

- Around the eyes and near the eyebrows (including temporal bone) – otherwise the sight may be permanently damaged!
- On the neck; the forehead up to the skulls and on the hairy head skin,
- On the nipples and areolas,
- On the genitals and their near surrounding areas,
- On the anus and its surrounding areas,
- Inside of the nostrils and ears,


- On particularly sensitive areas (such as ankles, bony areas etc.)
- On the areas you wish to have hair on in the future.

Moreover, men must not use the device on the face, neck and on the beard (as a variation of shaving machine or beard trimmer). In such case (the beard or facial hair) the epilation may result in lasting effects - uneven results, which may cause undesirable changes in characteristic facial features.

Hair removal using IPL technology may cause increased hair growth in some individuals. Based on the data available, the highest risk areas are the facial and neck ones, and it concerns rather women of Mediterranean, Middle Eastern and South Asian types.

Dear customers, please also note that regarding the specific principle of IPL technology operation and an individual character of human organism response, the seller cannot guarantee the full extent of the effects of the device (epilation, as well as rejuvenation functions). There may be a group of individuals, in which IPL treatment is less effective than in other. The period, within which the primary effect can be seen, may also be very individual – from approx. 4 weeks up to need to repeat the treatment in approx. 12-18 months.

Areas suitable for treatment

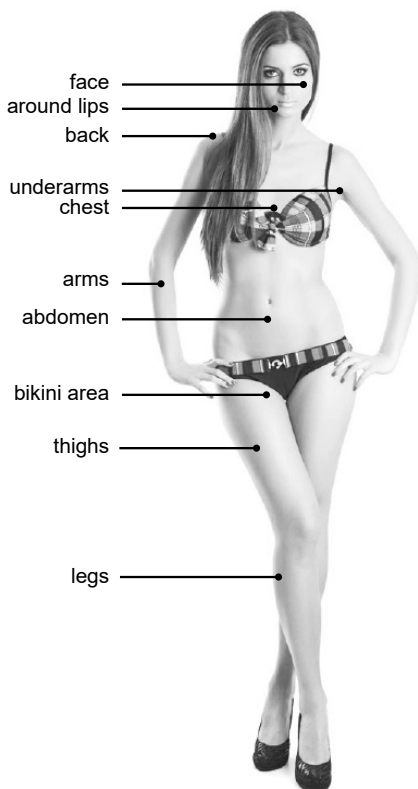III. FUNCTIONS OF DEVICE

The product is delivered with two attachments – yellow epilation attachment for hair removal and red attachment for skin rejuvenation (sometimes also photo-rejuvenation).

EPILATION

How does the IPL epilation work?

If used properly, the IPL technology (intensive pulses of light) epilation allows to achieve long-lasting effect and to prevent the regrowth of the hair. It is achieved by heating the hair and its root under the surface of the skin – the emitted light is absorbed, leading to destruction of their cells without damaging surrounding tissues. Consequently, the hair sheds within 1-2 weeks (Note: there may an impression of hairs continuing to grow, but they are only pushed out and then they fall out naturally) and its further regrowth on given area is prevented. If any new hairs occur (e.g. due to age or hormonal changes), they are fine and soft, and their density is reduced.

The major role plays the dark pigment (melanin), which can absorb the flashes of the light. The darker the hair is, the more amount of light it can absorb. That is why its use is not effective in light hair – as light hair absorbs only a minimum of light; majority of hair remains undamaged.

However, IPL epilation is only effective when the hairs are in their **growing phase!** There are approximately 10-30 % of them at the time of the session (this is very variable, and it depends on the individual, as well as the part of the body). Therefore, an effective hair removal requires some patience and repeating the epilation in regular cycles within a prolonged period (see below).

SKIN REJUVENATION

The flashes during the treatment generate heat inside the skin, stimulating the creation of collagen, which improves the elasticity and smoothness of the skin. At the same time, the heat treatment positively influences undesirable hormonal pigmentation (acquired, for example, during pregnancy), small veins, and smoothens various scars. The colour of the skin is unified to certain extent, its texture improved, the wrinkles are partially smoothed and overall freshness and vitality of the skin is improved.

For whom is the photorejuvenation suitable?

Due to hormonal changes, the attachment for rejuvenation is not recommended for use by persons under 25 years of age. However, always consider your overall health condition and respect the warnings given in these instructions. Also, do not use the attachment to unify the areas of pigmentation with other than hormonal cause (e.g. pigmented nevi). If you have any doubts concerning the use on your specific pigmented areas, always consult your physician or dermatologist first.

IV. DEVICE AND ACCESSORIES DESCRIPTION (picture 1)

A – IPL epilator

- A1 – ON./OFF button./flash intensity adjustment
- A2 – light activation button
- A3 – air vent grate
- A4 – air intake grate
- A5 – device socket
- A6 – light-signalling of scanning

B – LCD display

- B1 – marking of placed attachment
- B2 – set flash intensity
- B3 – remaining amount of flashed (longevity) of the attachment

C – IPL attachments

- C1 – hair removal attachment (HR); (order code ETA834700010)
- C2 – skin rejuvenation attachment (SR); (order code ETA834700020)
- C3 – application window with protection filter
- C4 – skin sensor

D – Safety eyewear (obj. kód ETA834700030)

E – Adapter

V. INSTRUCTION FOR USE

PREPARING FOR USE

To achieve optimal results, it is advisable to pre-treat the areas you wish to use the device on.

A) Pre-treatment of the skin on the area to be treated

- 1) Remove hair on the surface of the area in usual way (e.g. by shaving); if you have been recently epilated by a method including root removal, before using IPL, wait until the hair regrows. Using IPL technology on the areas, which are not hair-free, may result in undesired effects. e.g. burning of hair accompanied by smell, or contamination or damage to the application window.
- 2) Wash the skin, make sure it is free from any cosmetic products, perfumes etc. and let it dry.
- 3) Check, and, if needed, clean the application window with protection filter **C3**. It must not be contaminated with caught hairs, dust or any other impurities.

B) Setting the light intensity

The device provides 5 different pulsed light intensities. The IPL method is more effective in higher intensity settings. However, individual specifics of each person need to be considered in accordance with the rules/warnings given in the instructions. If you have darker skin tone, use lower intensity setting (see part "Skin type" above). As well, if you experience any discomfort or pain in any part of your body, reduce the intensity.

After switching on the device, you can change the pulsed light intensity by pressing **A1** button. The degree set is displayed by respective number of lines **B2** on the display **B**.

NOTE

Light impulse may cause a sensation of warmth or heat. In some cases, your skin may also slightly tingle. However, it should never be painful.

C) Setting

Based on the warnings in these instructions, evaluate, if your current health and skin condition allow the use of this device.

HOW TO USE**Skin tolerance test**

Before using the device for the first time (regardless the intensity degree) it is necessary to perform the tolerance test on a small skin area. It will help you to verify, if this type of technology is well tolerated by you and to discover any possible unwanted reaction of your skin. The test should be performed at least 24 before a complete treatment.

The tested area should be hair-free, clean and free from any cosmetic products. Then release 2-3 flashes close to each other at a low intensity of 1-2. Wait 24 hours and check your skin for any unwanted reactions (erythema, the skin must not be sensitive or swollen, nor you should experience any discomfort). If any complications occur, do not use the device and for further actions consult your physician/dermatologist.

After having performed the skin test, do not perform any treatment on the tested area for at least one week (by any of the attachments).

Modes of use

2 basic modes of the use of the device are available.

a) **manual** - it is a default mode, which can be simply described as a mode when each flash of light must be induced by pressing activation button **A2**. Use for more precise treatment of your face or other, rather small, parts of your body.

Proper contact with your skin and the possibility to release a flash is signalled by flashing **A6** signal lamp and **C4** sensor. If the application window is not in proper contact with your skin, the **A6** signal lamp and **C4** sensor are lit.

b) **automatic** - in this mode you can just press the application window **C3** onto your skin with no need to press the activation **A2** button. Then the flash releases in certain intervals will be activated automatically. It is suitable for treatment of larger areas.

Active mode is signalled by continuous flashing of **A6** signal lamp and **C4** sensor. Once the application window gets close to the skin, the flashes start to be released automatically.

After switching on the device, the automatic mode is activated by pressing and holding the A2 button for 2 seconds. The switch-over is signalled by an acoustic signal. Use the same procedure for switching the device back to the manual mode.

WARNING

- At lower levels of setting, the intervals between the flashes (or heating the lamp) are very short. If you operate in the automatic mode, you must continuously slide the light exit window, to prevent multiple flashes on the same area.
- When the device is switched off and then switched on again, it is activated in the mode in which it was switched off. However, if the device has been unplugged, when switched on again, it is activated in the default manual mode.
- For safety reasons, if the application window does not move, if in manual mode, the release of flashes will be automatically interrupted after 8 seconds.

Using the device

- 1) Select and by simple snapping place the selected attachment. If the device is switched on without an attachment placed on it, the HR and SR will flash on the display **B** signalling the missing attachment. To remove the attachment, use the same procedure (see the detail on the picture 1).
- 2) Put the connector of the adapter to the device socket **A5**. Then put the adapter in the wall socket. The sign "OFF" will be lit on the display and an audible signal will sound. Now the device is in a standby mode.
- 3) To switch on the device, press and long hold the button **A1**. An audible signal will sound, the integrated vent will be activated, and the display will be lit, indicating the marking of the attachment placed (HR or SR), the flash intensity set and numeric indication of the amount of remaining flashes for given attachment
- 4) If needed, change the flash intensity setting and the mode of use.
- 5) Put on the safety eyewear **D**.
- 6) Press the device onto the skin area to be treated. In automatic mode, the flashes will be released immediately. Thus, it is necessary to continuously slide the application window to the next area, so there is only a partial overlap and always move to the new area of the skin. In manual mode, press and release the application button **A2** to release a flash.
- 7) When you have finished the treatment, press and long hold the button **A1** to switch off the device and unplug the adapter from the wall socket. In case of inactivity lasting for more than 5 minutes, the device will be switched off automatically. Then let the device and the attachment cool down. Then perform routine maintenance and store the device in well-ventilated place out of reach of children and persons with reduced mental capabilities.

Skin care after use

Approximately 20 minutes after the treatment, you can apply moisturizers or regenerative lotions to the treated area. If you experience discomfort or skin irritation, you can calm the area by applying a cooling cloth (for 30-60 minutes). Do not use irritating cosmetic products! Do not moist the treated area with water for approximately 2 hours after the treatment.

Tips and advice on use

- Do not miss a place – release the flashes close to each other.
- When treating the face (women only!), more attention should be paid, as the skin in this area is very sensitive. Before application, close your eyes or look away (even with safety eyewear) and avoid releasing double flash on the same area. When treating the area over the upper lip, avoid the contact with lips and nostrils, which may be particularly sensitive.
- It is not recommended to use the device in a dark room. At a reduced light levels the flashes are more distinctive and may be exhausting to your eyes.
- During the treatment (but not during the session!) and for at least two weeks after the last treatment, protect your skin against sunlight using sunblocks SPF 30 and higher.
- Do not bath in disinfected pools etc. for at least 2 days.
- If a stronger skin reaction than a slight redness occurs, interrupt the use of the device and read "POSSIBLE SIDE EFFECTS".

TREATMENT SCHEDULE

Since the resulting effectiveness of the IPL epilator may vary between different individuals, please note that the length of successful treatment period may also differ. One will need only few treatments (lasting e.g. for 2 months), another will need a year or more. During the use, follow the treatment frequency as per following rules:

1st and 2nd month – use once a week.

3rd month – use 2 times a month (approx. after 15 days)

4th month – 1 year – use once a month

After 1 year – 1-2 times a year, if needed

VI. POSSIBLE SIDE EFFECTS

In principle, all methods of epilation and cosmetic interventions to a human body may cause some adverse effects. IPL technology is not an exclusion. These are mostly common reactions, which disappear within few hours/days. In case of more severe reactions (e.g. excessive redness or more intensive pain as a result of not following the rules in these instructions (mainly as a result of setting a high intensity, repeated application on the same area or application on a skin which is not hair-free) it is necessary to interrupt the treatment, evaluate the correctness of the procedure in accordance with the rules in the instructions and in case of any doubts or unusual skin reaction consult the further use of the device with your dermatologist or physician.

Common and harmless side effects:

- **Mild pain** – pain is perceived differently by different people. Under normal circumstances, the treatment with IPL device should not be painful. However, regarding different pain thresholds, it is possible that some people may experience certain discomfort (e.g. burning, prickle, tingling). But the pain should not be intensive and should disappear quickly. If you experience an excessive discomfort or pain at lower settings, do not use the device and consult your physician/dermatologist.

- **Slight redness** – this reaction is relatively common and in most cases it disappears quickly. It may occur immediately or within 24 hours after the treatment. The skin may show similar reaction as after sunburn. It is recommended to perform the treatment before sleep, so the possible redness does not interfere with your daily routine. However, if the problems persist for more than approx. 3 days, consult the further use of the device with your physician/dermatologist.
- **Smell of burnt hairs** – it usually occurs if you have not removed your hair before use and the flash burned the upper part of hairs above the skin surface. In such case stop using the device and remove hair from the treated areas.
- **Sensation of dryness and itching** – especially in combination with mechanical hair removal. It is common and it will disappear within few days. You can ease the symptom (but not immediately after use!) using a non-scented moisturizer.

Rare side effects:

- **burns, excessive redness or swelling** – often as a result of using a light intensity that too high for your skin tone.
- **discoloration** – it may manifest itself as either a darker patch or a lighter patch than the surrounding area. This also often occurs as a result of using the intensity that is too high for your skin tone.
- **Skin inflammation** – various infections and inflammations, that may occur after using the device on wounded areas (e.g. by shaving) or on ingrown hairs.

If the above-mentioned problems (excessive redness, discoloration, inflammation around hair follicle) do not disappear within 3 days, consult your doctor – some problems may be induced by other causes constituting a separate health problem. Do not use the device until your skin regenerates or restores its original colour.

VII. MAINTENANCE

To ensure an optimal condition of the device, clean it regularly after each session (or during the treatment, if needed). Neglecting the maintenance may lead to reduced effectiveness of the device.

Before performing maintenance, let the device and the attachment used cool down. Do not use abrasive, aggressive or flammable cleaning products! Wipe the surface of the device with a soft moist cloth. Check the vent slots and, if needed, remove e.g. caught hairs or other impurities. After each use, check the used attachment and its application window for any impurities caught on it. To clean the attachment and the application window, use only a dry soft cloth, not shedding fibre. The device itself and its attachments must not be washed or immersed in water or other liquids!

Storage

Put the appliance with the accessories back to the original package at a dry and dust free place, out of reach of children and incapacitated people. Ensure the protection of the application window with glass filter against damage.

VIII. ENVIRONMENTAL PROTECTION

If the dimensions allow, there are marks of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal, hand them over at the special collection places where they will be received free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal. Ask for more details at the local authorities or at a collection site. Fines can be afflicted for improper disposal of this type of waste in agreement with the national regulations. The appliance is equipped with rechargeable Li-ion batteries with long life. Dispose of the batteries in a suitable and safe way through the special collection networks. Never dispose of the batteries by burning!

More extensive maintenance or maintenance requiring intervention in the inner parts of the appliance must be carried out by professional service! Failure to comply with the instructions of the manufacturer will lead to expiration of the right to guarantee repair!

IX. TECHNICAL DATA

Voltage (V)	specified on the type label
Trimmer weight (kg)	0.2
Appliance Protection Class	III.
Adapter Protection Class	II.
Dimensions (h x l x w)	185 x 155 x 80 mm (with attachment)
Application area	approx 4 cm ²
Light spectrum	510-1100 nm
Noise level: Acoustic noise level of 37 dB(A) re 1pW.	

The manufacturer reserves the right to change the technical specifications and accessories for the respective models.

WARNINGS AND SYMBOLS USED ON THE APPLIANCE, PACKAGING OR IN THE INSTRUCTIONS MANUAL:

HOUSEHOLD USE ONLY. DO NOT IMMERSE IN WATER OR OTHER LIQUIDS.

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

 NOTICE

Do not use the appliance close to bath tubs, showers, bath basins or other water containers.

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© DATA 10/7/2019