

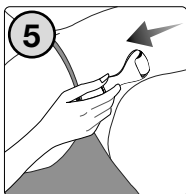
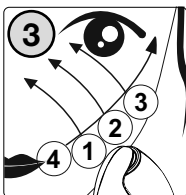
<i>Vibrační masážní váleček</i> • NÁVOD K OBSLUZE	CZ 4-7
<i>Vibračný masážny valček</i> • NÁVOD NA OBSLUHU	SK 8-11
<i>Vibrating massager roller</i> • USER MANUAL	EN 12-16
<i>Vibrációs masszázshenger</i> • KEZELÉSI ÚTMUTATÓ	HU 16-19
<i>Wibracyjny wałek masażowy</i> • INSTRUKCJA OBSŁUGI	PL 20-23
<i>Viebrierender Massageroller</i> • BEDIENUNGSANLEITUNG	DE 24-27

Fenité

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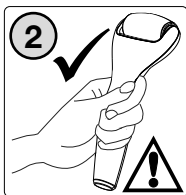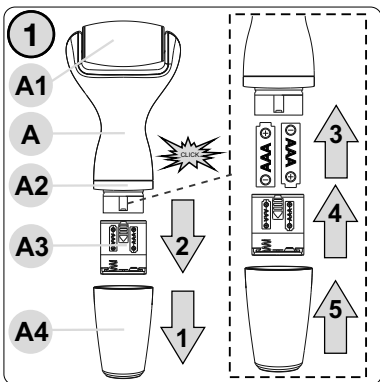Symbol  znamená UPOZORNENIE

Vibrating massager roller**eta 2353****INSTRUCTIONS FOR USE**

Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the warranty, the receipt and, if possible, the box with the internal packing.

I. SAFETY INSTRUCTIONS

- Consider the instructions for use as a part of the appliance and pass them on to any other user of the appliance.
- The appliance is only designated to be used in household and for similar purposes. It is not designed to be used in beauty parlours or for any other commercial utilisation!
- **Children and persons with limited physical or mental capabilities or lack of experience may only use this appliance if they are supervised or have been instructed of the use of the appliance in a safe manner and understand possible risks. Unattended children must not perform cleaning and maintenance intended to be performed by users. Children must not play with the appliance.**
- Never use the appliance if it does not correctly work or fell on the ground and is damaged.
- Before use make sure to remove glasses, contact lenses, and/or jewellery. We advise to tie your hair before any use. Make sure hair and loose pieces of clothing do not get into the roller to avoid risk of them tangling.
- When handling the appliance, please proceed so that no personal injury or damage to the appliance will occur.
- Never use the appliance if it does not correctly operate, fell on the ground or suffered any damage.

- Protect the appliance from dust, chemicals, excessive humidity, toiletries, aggressive liquid cosmetic preparations and intensive temperature changes, and do not expose it to direct sunlight.
- Do not insert the appliance in any bodily orifices.
- Never use the appliance for other than the intended purposes as described in these instructions for use!
- The manufacturer is not responsible for any loss caused by improper handling of the appliance and its accessories (e.g. injury, harm to health) and is not responsible according to the appliance warranty in the case of non-compliance with the safety instructions above. This appliance is not intended as a substitute for medical care.
- Use zinc-carbon or alkaline battery of AAA size. Observe the indicated polarity when exchanging batteries. If you are going not to use the Massager roller for a long time, remove the batteries.
- If the battery shows any indications of corrosion or leaking electrolyte, immediately replace it, otherwise the massager roller might get damaged.
- Dispose the empty battery in the appropriate manner (see article V. ECOLOGY).
- Keep the battery and Massager roller out of reach of children and incapacitated persons. Person who has swallowed the battery must immediately seek for medical assistance!
- **WARNING:** There is a risk of injury in case of incorrect use (not in accordance with the manual) of the appliance.
- The producer does not undertake any responsibility for damage caused by improper use of the appliance (e.g. injury etc.) and its guarantee does not cover the appliance in the case of failure to comply with the safety warning above.

SPECIAL SAFETY INSTRUCTIONS

- The device must not be used by: Pregnant or breastfeeding women and If your skin is very sensitive.
- Make sure to use the appliance only on face and body skin, not on eyes.
- Do not use the appliance:
 - if your skin is oversensitive, damaged or exematic (e.g. in the case of fresh scars, small injuries, burns, erythema, swollen cheeks, skin diseases...). To avoid dissemination of micro inflammations, do not use the massager roller including cases when your skin suffers acne. We do not recommend using the massager roller on injection-treated skin provided that in the given area you have undergone surgical intervention or plastic facial rejuvenation (such as e.g. Botox, facial fillers, thread lifts etc.).
 - On areas which may present a health hazard (varicose veins, moles, birthmarks).
 - If you have any tattoo on your face or permanent makeup.

- If you feel pain or experience any unpleasant feeling during the using, or if any other health problem appears, immediately switch off the appliance, and consult its further use with your practitioner or dermatologist.
- The appliance is not intended to cure any skin disease symptoms. Please consult use of the appliance with your doctor, where necessary.
- Do not use the appliance on children, incapacitated persons or persons unable to express their feelings.

Note

ETA a.s. is not a licensed provider of health and cosmetic care and does not have expertise in the field of treatment, diagnostics of the health condition and related effects or suitability for the use of appliances that may affect the health condition in any way. During any exercise, massage therapy or cosmetic intervention, there is a theoretical possibility of physical injury or, for example, worsening of the symptoms of existing diseases. By using this appliance, you agree that you use it at your own risk and assume all risks of injury or deterioration. You also agree to waive any claims (whether known or unknown) against ETA a.s. arising from improper use of this appliance or through negligence. Before using any appliances that may affect your health, first consult your doctor about the suitability of their use.

II. DESCRIPTION OF FEATURES (PIC. 1)

A – Massager roller body

A1 – roller

A3 – battery compartment

A2 – sensor

A4 – cover

III. INSTRUCTIONS FOR USE

Remove all packaging materials and take out the appliance.

INSERTING THE BATTERY

Remove the cover **A4** and battery cover **A3** as per Fig. 1 and insert a new AAA battery (2 pcs.) into the compartment **A**. Close the cover in the opposite way.

Use

Hold the appliance body and make sure that your fingers are in full contact with the sensor (Pic. 2). Move the massage roller over the skin, the appliance will start to vibrate automatically. If you move the appliance away from the skin, the vibrations will be switched off. Face: gently move the massage roller from the jaw to the cheekbone (Pic. 3). Then slowly move the roller from the mouth to the ear (Pic. 3). Neck: move the massage roller from the shoulder along the line of the neck and jaw (Pic. 4). Arms: Gently move the roller from the wrist to the collarbone (Pic. 5).

IV. MAINTENANCE

For maximum hygiene, clean the roller **A1** after each use (wipe the massage head surface with a wet soft cloth)! Do not use any abrasive or aggressive cleaning agents! Do not immerse the massager body **A** in water. Store the appliance at a safe, dry and dust-free place, out of reach of children and incapacitated persons.

V. ENVIRONMENT

If dimensions permit, all pieces have printed symbols of materials used for production of packing, components and accessories as well as recycling symbols. The symbols on the product or documentation means the electric or electronic parts used in the appliance may not be disposed off together with communal waste. Hand over the electric or electronic part to a waste collection yard for free disposal. Proper disposal of the product will help preserving valuable natural resources and preventing potential negative impacts on environment and human health that would otherwise result from improper waste disposal. For more details, please contact your local authorities or collection yard. Fines may be imposed for improper disposal of this type of waste pursuant to national directives. Remove battery from the appliance and dispose it off in a suitable way using special collection centres. Do not incinerate the battery!

VI. TECHNICAL DATE

Protection class of the appliance III.

The manufacturer reserves the right to change the technical specifications for the respective models. Expert service must do the large maintenance or to intervene into inside parts of the appliance! Non-observance of the instructions of the manufacturer cancels the right for guarantee repair!

WARNINGS AND SYMBOLS USED ON THE APPLIANCE, PACKAGING OR IN THE INSTRUCTIONS MANUAL:

HOUSEHOLD USE ONLY.

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO

NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

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© DATE 20/5/2022