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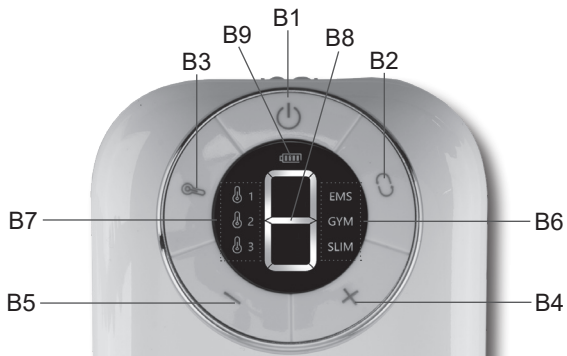
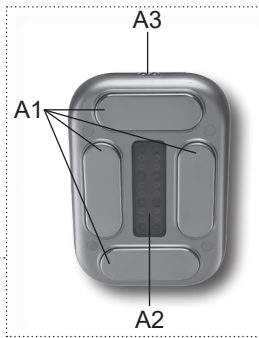
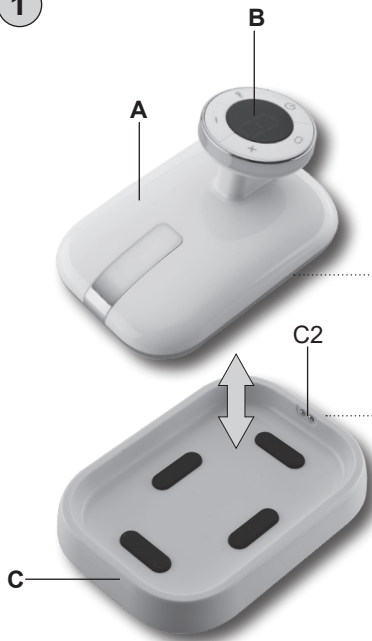
**HOUSEHOLD USE ONLY. DO NOT IMMERSE IN WATER OR OTHER LIQUIDS.**

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Fenité

## INSTRUCTIONS FOR USE

Dear customer, thank you for purchasing our product. Keep these instructions including the warranty, the receipt and, if possible, the box with the internal packing.

## I. SAFETY WARNING

- Please consider the instructions herein as a part of the appliance and pass them over to the next user of the appliance.
- **The appliance is only designated to be used in household and for similar purposes. It is not designed to be used in beauty parlours or for any other commercial utilisation!**
- This appliance may only be used by children aged 8 years or older, and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- **Never immerse the appliance in water!**
- **When handling the appliance, please proceed so that no personal injury or damage to the appliance will occur.**
- Protect the appliance from dust, chemicals, excessive humidity, intensive temperature changes, and do not expose it to direct sunlight.
- Before use, take off all metal items and jewellery (e.g., glasses, watch, earrings, bracelets, rings, etc.).
- Do not insert the appliance or its accessories in any bodily orifices.
- Never use the appliance for other than the intended purposes as described in these instructions for use!
- Only use the appliance with the accessories designed for this type. Use of unapproved accessories may cause injury.
- This appliance contains batteries that may only be replaced by trained personnel.
- This appliance has a heated surface, and the persons with impaired heat perception must be careful in using it.
- **WARNING:** In the case of incorrect use of the appliance that contradicts the operating instructions, there is a risk of injury.
- Texts in a foreign language and images provided on the packaging or on the product, if any, are translated and explained at the end of this language version.

- The manufacturer is not responsible for any loss caused by improper handling of the appliance and its accessories (e.g. injury, harm to health) and is not responsible according to the appliance warranty in the case of non-compliance with the safety instructions above. This appliance is not intended as a substitute for medical care.
- Do not use the appliance if your skin is oversensitive, damaged or exematic (e.g. in the case of fresh scars, small injuries, burns, erythema, swollen cheeks, skin diseases...).
- If you perceive pain/ unpleasant feelings or if any other health problems occur, switch off the appliance immediately and consult with your physician or a dermatologist about its use.
- The appliance is not intended to cure any skin disease symptoms. Please consult use of the appliance with your doctor, where necessary.
- Do not use the appliance on children, incapacitated persons or persons unable to express their feelings.

### CONTRAINDICATIONS

The device **MUST NOT** be used by persons affected by any of the following contraindications: diabetes, cancer, epilepsy, viral diseases, high fever, liver disease, kidney disease, skin disease, high cholesterol, fatigue syndrome, heart disease, pacemaker, metal implants in the body, vascular disorders, pregnancy, use of anti-clotting drugs, organ transplants.

### FORBIDDEN AREAS

The device **MUST NOT** be used in the following areas: face area, heart area, thyroid area, armpit, elbow, wrist, hand, knee, ankle, foot, spine area.

## SPECIAL SAFETY INSTRUCTIONS

### WARNING

- To power the appliance, only use a source of safe low voltage of 5V which complies with applicable standards (mobile phone chargers with a micro USB connector with double insulation are suitable).
  - Do not use USB ports on other devices (e.g. computers, laptops, tablets, mobile phones, etc.) to power the appliance. Such devices may not be dimensioned for the current consumption of the appliance and it could damage them.
- Children must be kept away from the appliance and its adapter.
  - The appliance has no charging source of its own, it is only designed for charging via USB (an adapter with USB output with output voltage of 5.0 V, and output current of at least 1000 mA). Before making any connection, please pay attention to instructions of the manufacturer of the respective power source.

- Never use the appliance if the adapter is damaged, if the appliance does not operate properly, fell on the ground and suffered any damage, or fell into water. In such case bring the appliance to the specialised service centre to verify the appliance safety and proper operation.
- For additional protection we recommend installing a protective switch (RCD) in the electrical circuitry of the bathroom, with the nominal actuating current not exceeding 30 mA. Ask an inspecting technician or an electrician for advice.
- The appliance and the connected adapter must not be used at places where it could fall into bath tub, wash basin or swimming pool. If the appliance fell into water anyway, do not take it out! First of all, unplug the adapter from the socket and only then take the appliance out. Then bring the appliance to a specialised service centre to check its safety and proper functionality.
- If the appliance is used in a bathroom, it must be disconnected from the mains after use by unplugging the power cord from the wall socket, because the appliance poses a risk if near to water, even when switched off.
- Always unplug the appliance from charger if you leave it unattended and cleaning.
- The mains adapter supplied with the appliance is only designated for this appliance.
- Never plug the adapter in the electrical socket or unplug it with wet hands.
- Never immerse the adapter into water or other liquids (not even partly) and do not wash it in running water!
- The appliance gets hot when charging, which is absolutely normal.
- Charge the battery at a temperature close to the room temperature.
- Never charge the appliance in direct sunlight or close to heat sources, and if the temperature is below 10 °C or above 40 °C, adversely affects the accumulator service life.
- Prevent hanging the power adapter cord loosely over desktop edge where children may reach it.
- The power cord of the adapter must not be damaged with sharp or hot objects or open fire, and must not be immersed into water or bent over sharp edges.

**Note**

ETA a.s. is not a licensed provider of health and cosmetic care and does not have expertise in the field of treatment, diagnostics of the health condition and related effects or suitability for the use of appliances that may affect the health condition in any way. During any exercise, massage therapy or cosmetic intervention, there is a theoretical possibility of physical injury or, for example, worsening of the symptoms of existing diseases. By using this appliance, you agree that you use it at your own risk and assume all risks of injury or deterioration. You also agree to waive any claims (whether known or unknown) against ETA a.s. arising from improper use of this appliance or through negligence. Before using any appliances that may affect your health, first consult your doctor about the suitability of their use.

**II. DESCRIPTION OF THE APPLIANCE (Fig. 1)****A – Appliance body**

- A1 – application areas
- A2 – red LED light for photon therapy
- A3 – charging contacts

**B – Control panel with LED display**

- B1 – on/off button + vibration adjustment button
- B2 – mode setting button
- B3 – button for setting the heating level
- B4 – "+" button to increase intensity
- B5 – "-" button to reduce intensity
- B6 – set mode indicator
- B7 – set temperature indicator
- B8 – display of the set intensity level
- B9 – battery level indicator

**C – Recharging base**

- C1 – charging port
- C2 – charging contacts

**D – USB-C cable \*****E – fixation belt on the arm****F – fixation belt on the body**

\* USB-C cable: A source with USB-C terminal may be connected to USB cable (for example AC power adapter, backup power source, etc.).

### III. INSTRUCTIONS FOR USE

Before using the appliance, its battery must be charged first. Charging is possible with the enclosed USB cable from the mains, i.e. an AC adapter, or a backup power source.

#### CHARGING

- Before the first use or if the appliance was not used for a long time, charge it completely.
  - The standard charging time is approximately 3 hours. and the fully charged accumulator allows for the operating time of up to 12 hours. The operating time decreases depending on the appliance loading and due to the accumulator ageing and wear and tear.
1. Place the appliance on the charging base **C** as shown in Fig. 1 and connect it to the power supply (charger) **D** via the USB-C cable. The appliance cannot be used while charging.
  2. Do not exceed the designated charging time!

#### Note

Battery level and charging state are displayed by indicator **B9** (flashing = charging; steady light = fully charged battery).

#### CONTROL

##### Control panel with LED display (B)

Button		Description of function
Marking	Symbol	
B1		Long pressing - switching the appliance on/off; Short pressing - on/off vibration and switching between different vibration levels (3 levels - low, medium, high).
B2		A short press switches between the different modes of operation (3 modes - EMS, GYM, SLIM). The set mode is displayed on the <b>B6</b> indicator.
B3		A short press switches between the different heating levels (3 levels - low, medium, high). The set temperature is displayed on the <b>B7</b> indicator.
B4	+	Briefly press the corresponding button to increase/decrease the intensity level (9 intensity levels). The set intensity is displayed on the <b>B8</b> indicator.
B5	–	

### Comment

Each press of the button is accompanied by an acoustic signal.

## USE OF SONOGEL

In order for the appliance to work properly, it is necessary to apply a water-based sonogel to the area (parts) to be treated. This gel forms a necessary bonding layer between the application surfaces and the skin. This allows the EMS impulses to better penetrate the skin. Sonogel is not included in the package but is available at any pharmacy or medical supply store.

## DESCRIPTION OF FUNCTIONS AND THEIR USE

- Regular use of the appliance helps to break down fatty tissue, stimulates muscles, increases skin elasticity and counteracts cellulite.
- The appliance is particularly suitable for the treatment of arms, thighs, calves, tummy and buttocks.

### Vibration

They gently massage to stimulate the skin and give a pleasant sensation during application.

### Red light

Red light photon therapy improves blood circulation in the skin, promotes natural collagen production and reduces pigment formation.

### EMS mode

It uses a micro current to stimulate muscle contractions.

### GYM and SLIM modes

Both modes use electrical impulses to stimulate the muscles (contraction + relaxation).

## APPLICATION

- First apply a layer of sonogel to the area (parts) to be treated.
- Switch on the appliance and select the desired operating mode.
- Apply the appliance to the area with the **A1** application surfaces and pass over the area with slow movements and gentle pressure. **E** and **F** belts can be used for fixation (on the arm or in the tummy area).
- Switch off the appliance when you have finished working.

### WARNING

Use the appliance for a maximum of 10 minutes and once a day in each part.

## AUTOMATIC SWITCH-OFF

After about 10 minutes of operation, the appliance always switches off automatically.

## IV. MAINTENANCE

Only wipe the appliance, charging base, fixation belts and USB-C cable with a soft, damp cloth. Do not use any abrasive or aggressive cleaning agents! Guarantee of 24 months does not apply to decrease in the accumulator's effective capacity caused by its use or age. Effective capacity gradually decreases depending on the method of use of the accumulator.

**Expert service must do the large maintenance or to intervene into inside parts of the appliance! Non-observance of the instructions of the manufacturer cancels the right for guarantee repair!**

## V. ENVIRONMENT

If the dimensions allow, there are labels of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal, hand them over at special collection places where they will be accepted free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal. Ask for more details at the local authorities or at a collection site. Fines can be imposed for improper disposal of this type of waste in agreement with the national regulations. Dispose of the battery safely by means of collection networks intended for this purpose. Never dispose of the batteries by burning!

### Taking out the battery

Only take out the battery from the appliance after it is fully discharged, and after the end of the appliance's service life! In the course of operation, the battery does not need to be removed! For taking out the battery, remove rubber screw covers using a screwdriver, unscrew the screws and take out the battery.

## VI. TECHNICAL DATA

Voltage (V)	shown on the type label of the appliance
Absorbed power (W)	shown on the type label of the appliance
Weight (kg) approximately	0,2
Protection class of the appliance	III.
Size, (mm)	90 x 131 x 71

The manufacturer reserves the right to change the technical specifications and accessories for the respective models.

### WARNINGS AND SYMBOLS USED ON THE APPLIANCE, PACKAGING OR IN THE INSTRUCTIONS MANUAL:

HOUSEHOLD USE ONLY. DO NOT IMMERSE IN WATER OR OTHER LIQUIDS.

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, ARRANGES OR PLAYPENS. THIS BAG IS NOT A TOY.

WARNING

Do not use the appliance close to bath tubs, showers, wash basins or other water containers.

Please read the instructions manual

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© DATE 31/05/2023