

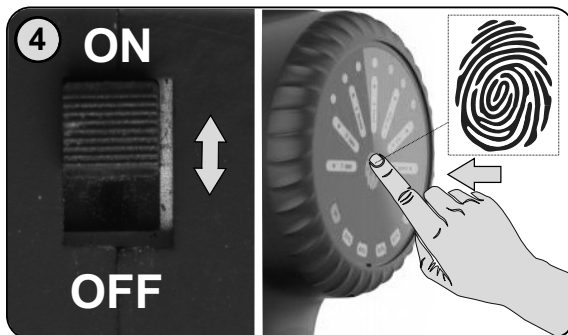
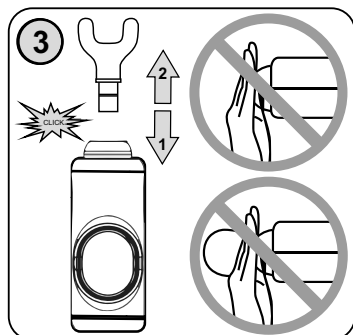
<i>Masážní pistole</i> • NÁVOD K OBSLUZE	CZ	3-7
<i>Masážna pištoľ</i> • NÁVOD NA OBSLUHU	SK	8-13
<i>Massage gun</i> • USER MANUAL	EN	14-18
<i>Massázspisztoly</i> • KEZELÉSI ÚTMUTATÓ	HU	19-23
<i>Pistolet masażowy</i> • INSTRUKCJA OBSŁUGI	PL	24-29
<i>Massagepistole</i> • BEDIENUNGSANLEITUNG	DE	30-35

Power Care Pro

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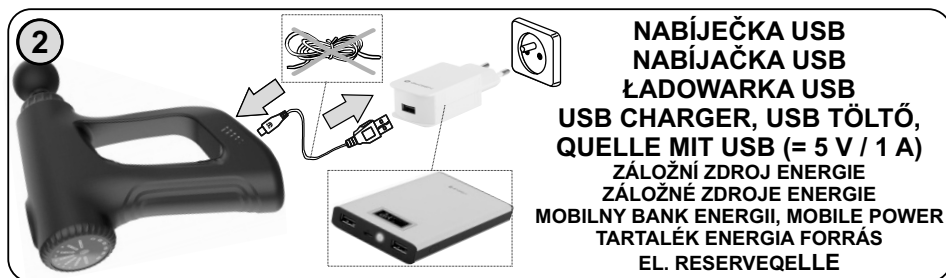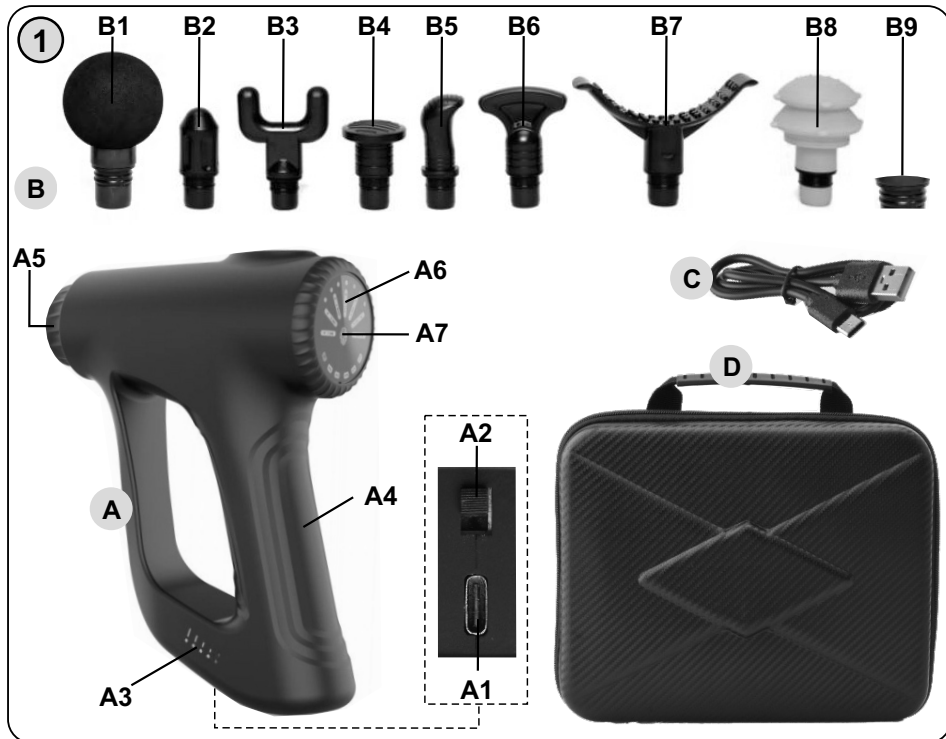HOUSEHOLD USE ONLY | Len na použitie v domácnosti.

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, ARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

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Power Care Pro

INSTRUCTIONS FOR USE

Dear customer, thank you for purchasing our product. Please keep these instructions including the receipt and, if possible, the carton box with the internal packing.

I. SAFETY WARNING

- Please consider the instructions herein as a part of the appliance and pass them over to the next user of the appliance.
- The appliance is only designated to be used in household and for similar purposes. It is not designed to be used in massage saloon or for any other commercial utilisation! The appliance is not a medical instrument and is not used for medical purposes.
- This appliance may only be used by children aged 8 years or older, and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Never immerse the appliance in water!
- When handling the appliance, please proceed so that no personal injury or damage to the appliance will occur.
- Protect the appliance from dust, chemicals, excessive humidity, toiletries, agresive liquid cosmetic **preparations and intensive temperature changes, and do not expose it to direct sunlight.**
- Do not use the appliance during bathing or showering and do not apply it to wet skin. The Massage gun is intended for dry use only.
- Do not insert the appliance or its accessories in any bodily orifices.
- Never use the appliance for other than the intended purposes as described in these instructions for use!
- Only use the appliance with the accessories designed for this type. Use of unapproved accessories may cause injury.
- **WARNING:** In the case of incorrect use of the appliance that contradicts the operating instructions, there is a risk of injury.
- Texts in a foreign language and images provided on the packaging or on the product, if any, are translated and explained at the end of this language version.
- The manufacturer is not responsible for any loss caused by improper handling of the appliance and its accessories (e.g. injury, harm to health) in the case of non-compliance with the safety instructions above. This appliance is not intended as a substitute for medical care.

SPECIAL SAFETY INSTRUCTIONS

- Children must be kept away from the appliance and its adapter.
- The appliance has no charging source of its own, it is only designed for charging via USB (an adapter with USB output with output voltage of 5.0 V, and output current of at least 1 A). Before making any connection, please pay attention to instructions of the manufacturer of the respective power source.

- Never use the appliance if the adapter is damaged, if the appliance does not operate properly, fell on the ground and suffered any damage, or fell into water. In such case bring the appliance to the specialised service centre to verify the appliance safety and proper operation.
- For additional protection we recommend installing a protective switch (RCD) in the electrical circuitry of the bathroom, with the nominal actuating current not exceeding 30 mA. Ask an inspecting technician or an electrician for advice.
- The appliance and the connected adapter must not be used at places where it could fall into bath tub, wash basin or swimming pool. If the appliance fell into water anyway, do not take it out! First of all, unplug the adapter from the socket and only then take the appliance out. Then bring the appliance to a specialised service centre to check its safety and proper functionality.
- If the appliance is used in a bathroom, it must be disconnected from the mains after use by unplugging the power cord from the wall socket, because the appliance poses a risk if near to water, even when switched off.
- Always unplug the appliance from charger if you leave it unattended and cleaning.

WARNING

- **To power the appliance, only use a source of safe low voltage of 5 V which complies with applicable standards (mobile phone chargers with a USB-C connector with double insulation are suitable).**
- **Do not use USB ports on other devices (e.g. computers, laptops, tablets, mobile phones, etc.) to power the appliance. Such devices may not be dimensioned for the current consumption of the appliance and it could damage them.**

- The mains adapter supplied with the appliance is only designated for this appliance.
- Never plug the adapter in the electrical socket or unplug it with wet hands.
- Never immerse the adapter into water or other liquids (not even partly) and do not wash it in running water!
- The appliance gets hot when charging, which is absolutely normal.
- Charge the battery at a temperature close to the room temperature.
- Never charge the appliance in direct sunlight or close to heat sources, and if the temperature is below 10 °C or above 40 °C, adversely affects the accumulator service life.
- Prevent hanging the power adapter cord loosely over desktop edge where children may reach it.
- The power cord of the adapter must not be damaged with sharp or hot objects or open fire, and must not be immersed into water or bent over sharp edges.
- Do not use the appliance if you feel tired or weakened. After using painkillers, sedatives, narcotics or consumption of alcohol, the body sensitivity to pain is reduced.
- Do not use the appliance if your skin is oversensitive, damaged or exematic (e.g. in the case of **fresh scars, small injuries, burns, erythema, skin diseases...**).

- If you experience pain or discomfort during or after the massage, or if you experience any other health problems, switch off the appliance immediately and have a consultation with your doctor. Any bruises may occur regardless of the intensity setting or pressure applied. Check the frequently massaged parts and stop the massage immediately at the first signs of pain or discomfort.
- Talk to your doctor about the appropriateness of using the appliance if you are pregnant, have diabetes, neuropathy, heart failure, epilepsy, are injured or have undergone any type of surgery.
- Do not use the appliance on children, incapacitated persons or persons unable to express their feelings.
- Do not use the appliance in close proximity to the eyes, directly on the bones, head, spine or other sensitive parts of the body. Massage only soft (muscle) tissues.
- Keep your fingers, hair or other part of your body away from the shaft between the inlet hole and the inserted attachment, as this may pinch the skin / hair and genitals.
- Do not use the appliance continuously for more than one hour. Then let it cool for at least 30 minutes before using it again.
- This appliance contains batteries that may only be replaced by trained personnel.
- **The capacity of the product battery decreases over time. Decreasing of battery capacity occurs during normal use of the appliance and is considered usual wear and tear of the appliance. Decreasing battery capacity is not a defect of the appliance.**
- The accumulator must not be removed when the appliance is normally used. Only remove the accumulator if the appliance is to be discarded (see article V. ECOLOGY).

Note

ETA a.s. is not a licensed provider of health and cosmetic care and does not have expertise in the field of treatment, diagnostics of the health condition and related effects or suitability for the use of appliances that may affect the health condition in any way.

During any exercise, massage therapy or cosmetic intervention, there is a theoretical possibility of physical injury or, for example, worsening of the symptoms of existing diseases. By using this appliance, you agree that you use it at your own risk and assume all risks of injury or deterioration. You also agree to waive any claims (whether known or unknown) against ETA a.s. arising from improper use of this appliance or through negligence. Before using any appliances that may affect your health, first consult your doctor about the suitability of their use.

II. DESCRIPTION OF THE APPLIANCE (Fig. 1)

A – Massage gun

- | | |
|-----------------------------------|---------------------------------------|
| A1 – charging socket | A5 – inlet for extensions |
| A2 – ON / OFF button | A6 – display |
| A3 – light charging signalization | A7 – touch sensor intensity switching |
| A4 – handle | |

B – Massage attachments

- | | |
|------------------------|------------------------------|
| B1 – round extension | B5 – thumb attachment |
| B2 – conical extension | B6 – narrow wedge attachment |
| B3 – fork attachment | B7 – curved attachment |
| B4 – flat attachment | B8 – flat / air attachment |
| | B9 – distance ring (4 psc.) |

C – USB-C cable *

- * USB-C cable: A source with USB terminal may be connected to USB-C cable (for example AC power adapter, backup power source, etc.).

D – Protective case

III. INSTRUCTIONS FOR USE

Before using the appliance, its battery must be charged first. Charging is possible with the enclosed USB-C cable from the mains, i.e. an AC adapter, or a backup power source.

CHARGING (Fig. 2)

- Before the first use or if the appliance was not used for a long time, charge it completely.
- The standard charging time is **approximately 6-7 hours**, and the fully charged accumulator allows for the operating time of **up to 240 minutes**. The operating time decreases depending on the appliance loading and due to the accumulator ageing and wear and tear.

1. Make sure the appliance is switched off.
2. Connect the enclosed USB-C cable to the power supply terminal **A1**, and connect the other side to a power source, for example to a compatible AC power adapter, or to a backup power source. While charging, the appliance cannot be in use.
If the appliance is turned on and you start charging, the appliance shall switch off.
3. Do not exceed the designated charging time!
4. The status and progress of charging is indicated by display **A6** and by indicator **A3**.
5. Charge the battery as soon as you notice reduced power or the appliance switches itself off.

Operation signal lights and charging lights (A3, A6)

10x flashing / on continuously red light segment 	Battery discharged (the device will switch off)
LEDs A3 go off	Battery discharged (the device will switch off)
Increasing number of green flashing LEDs A3	Accumulator is charging.
Blue light segments A6 on continuously 	Accumulator is charged / indication of operation (individual segments turn off during operation, as the battery discharges)

Optimise the battery service life

- Do not leave USB-C cable connected to the appliance.
- Let the accumulator cool down prior to charging. During charging or regular operation, the appliance heats up, which is a normal phenomenon, and not a reason for a complaint.
- To extend the lifespan of the accumulator, we recommend discharging the accumulator completely 2x a year by normal use and then charging it fully again.

CONTROL (Fig. 4)

- Press the switch **A2** to the position „ON“ (indicator **B5** indicating connection to the mains). The massage gun will automatically switch to standby mode.
- You can switch between 10 vibration intensities by briefly pressing the **A7** touch sensor. The set intensity level is indicated by the corresponding indicator light on the display **A6**. When switched on, the first indicator light illuminates, indicating the lowest vibration intensity.
- By long pushing of button **A7** (approx. for 2 seconds) you switch the appliance on/off.

Notes: The massage gun is set to 1 intensity level when switched on 

USE

Choose the appropriate massage head to suit your needs and health to maximize the benefits and minimize the risks.

- You attach the attachments by simply sliding them (you can turn them slightly) into the inlet opening for the attachments **A5** (Fig. 3). You will then remove them in the opposite way.
- Only use the appliance on clean and dry skin (recommended through clothing). Do not use in combination with massage oils, etc.
- Place on the massaged area and gradually move the gun over the muscular areas with slight pressure. In general, you should not massage the same area for too long - preferably about 60 seconds for each area.

Type	Extension	Use
B1	round extensio	Universal gentle massage of larger areas, such as hands, calves, back, hips.
B2	conical extension	Targeted intensive massage of specific places such as palms, joints, feet.
B3	fork	Intensive massage of the neck, around the spine or eg Achilles tendons.
B4	flat	Universal intensive massage designed for most muscle groups and body parts.
B5	thumb	Suitable for smaller muscle parts for better focus and muscle relaxation. Painful areas of the feet, shoulders.
B6	narrow wedge	Suitable for painful spots in the lumbar/back area.
B7	curved	Suitable for shoulders and legs. Simulated tapping massage with constant intensity e.g. on arms, thighs and other curved parts of the body.
B8	flat / air	Suitable for „bones“ of different parts of the body. Larger surface area, the best to start with before concentrating the focus on the pain location.

Auto-off mode

After approximately 10 minutes of continuous operation the appliance automatically switches off.


IV. MAINTENANCE

The appliance must be cleaned after every use! Clean the appliance surface and the USB cable only with use of a soft damp cloth. First, take down the removable heads (**B1-B8**) of the appliance. Rinse the **B2-B8** attachments under running water and then wipe dry. Clean the foam attachment **B1** with a damp cloth, or you can rinse it (then dry it immediately). Do not use any abrasive or aggressive cleaning agents! If the distance ring **B9** on any of the attachments is worn out or damaged, replace it with a new one (the same type). Proceed in reverse order to assemble all the components. Store the massage gun at a safe and dry place, out of reach of children and incapacitated people.

More extensive maintenance or maintenance requiring intervention in the inner parts of the appliance must be carried out by professional service!

Failure to comply with the instructions of the manufacturer will lead to expiration of the right to repair!

The information regarding warranty and product repairs you can find on the website www.eta.cz/servis-eu.

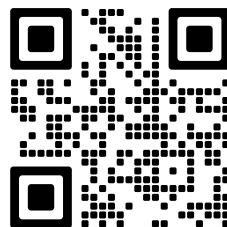V. ENVIRONMENT

If the dimensions allow, there are labels of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal, hand them over at special collection places where they will be accepted free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal. Ask for more details at the local authorities or at a collection site. Fines can be imposed for improper disposal of this type of waste in agreement with the national regulations. Dispose of the battery safely by means of collection networks intended for this purpose. Never dispose of the batteries by burning!

Taking out the battery

Only take out the battery from the appliance after it is fully discharged (i.e. let the appliance in operation until its motor stops), and after the end of the appliance's service life! The battery should not be removed and is not accessible during normal operation! With use of a suitable tool demount the appliance body and replace the accumulator. However pay attention to prevent any cut injury.

VI. TECHNICAL DATA

Voltage (V)	Specified on the rating place
Appliance weight (g)	779
Appliance protection class	III.
Dimensions (H x D x W) in mm	172 x 225 x 67
Noise level: Acoustic noise level of dB(A) re 1pW	45 - 60

The manufacturer reserves the right to change the technical specifications and accessories for the respective models.

Warnings and symbols used on the appliance, packaging or in the instructions manual

HOUSEHOLD USE ONLY, DO NOT IMMERSE IN WATER OR OTHER LIQUIDS

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, ARRANGES OR PLAYPENS. THIS BAG IS NOT A TOY.

It is not allowed to dispose it together with communal waste.

NOTICE

Please read the instructions manual

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© DATE 15/12/2023